

enovis™

# PRIMA SYSTEM

a LimaCorporate family product

## Efficient by design



**EFFICIENCY**  
Doing more, with less.

**TECHNOLOGY**  
Conceived within  
technology.

**PERFORMANCE**  
Built on a strong  
clinical heritage.

## PRIMA SYSTEM

# EFFICIENT BY DESIGN


**PRIMA TT Glenoid** features a new philosophy of glenoid replacement, enabling streamlined implantation with small monoblock implants for primary cases as well as providing modular wedged baseplates for defect-filling needs.

### 1. Tailored fixation

combining dia. 5mm peripheral screws with optional locking caps and central compressive screws, available in two diameters - 5mm and 6,5mm.

### 2. Streamlined glenoid implants

monoblock dia. 25mm baseplates with built-in TT central peg, featuring a one-step glenoid preparation. Regular and 10° Full Wedge options available.

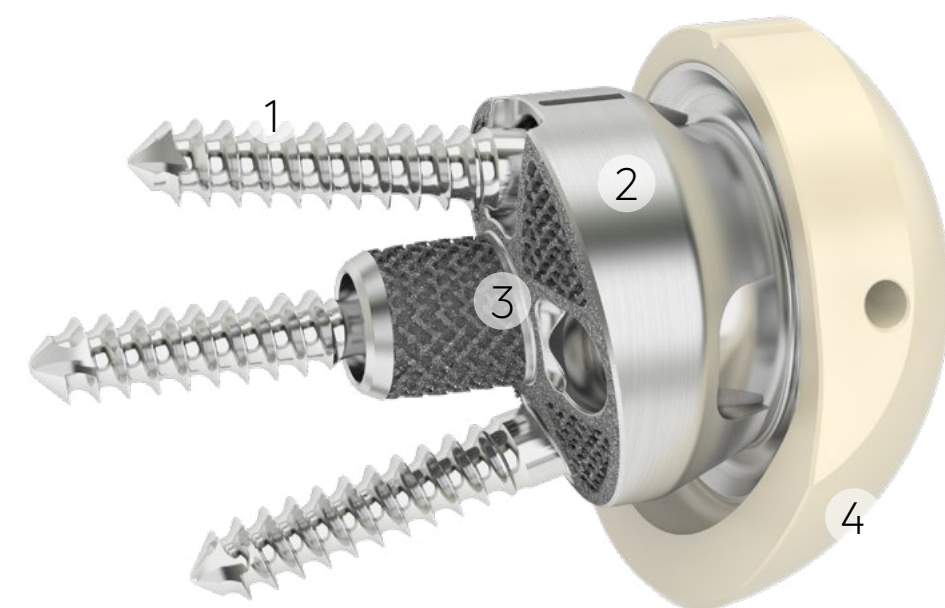

### 3. Defect-filling baseplates

modular Dia. 25mm and 28mm baseplates, featuring 4 different lengths of TT central peg to tailor glenoid fixation. Regular and 10°/15°/20° Full Wedge options available.

### 4. LimaVit glenospheres

unique design, available with two diameters: 40mm and 44mm. Both sizes feature low, medium and high CoR lateralization options.

**PRIMA Stem** is designed to address the need for efficiency and reproducibility. This 3D-printed bone-preserving solution enables the treatment of primary cases with a streamlined surgical technique and a reduced inventory.\*\*

**5. Inlay design** to allow soft tissue balancing and reduce risk of joint overstuffing.

**6. Proven Ring FiTT** 3D-printing technology and design to enhance metaphyseal fixation and osteointegration<sup>[2-4]</sup> where the best quality of bone is.<sup>[1,5-8]</sup>

**7. 5 sizes** of stem proximal outer diameter designed through FIGURA system to fit in a wide range of anatomies\*


**8. Fins** for rotational stability and **open design** to reduce the volumetric filling ratio<sup>[9]</sup> and to promote the load transfer on the **Ring FiTT**.

**9. Monoblock stem** with an anatomically-inspired 135° neck-shaft angle.



\* Data on file at Enovis. Simulation study does not necessarily indicate clinical performance.

\*\* Only one instrument set needed to implant PRIMA Stem in reverse configuration - humeral side only.

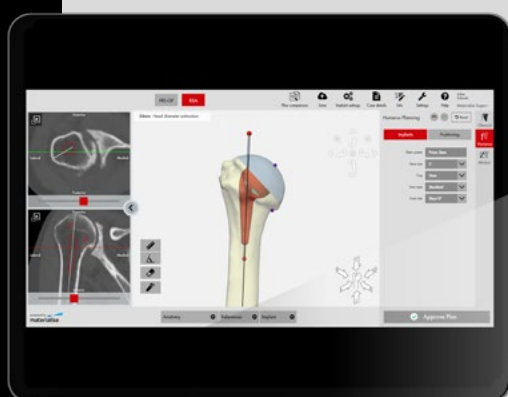

Find **PRIMA System** integrated in our Pre-Operative Planner software powered by Materialise\* to plan cases in a simple, automated and independent manner.

- Auto-Segmentation.
- User-friendly interface.
- 360° planning.
- Plan everywhere.
- Patient specific instruments.

\* Materialise NV is the proprietor and - to all effects - the legal manufacturer of the pre-op planning software.



[BIBLIOGRAPHY](#)

This publication is not intended for distribution in the U.S.  
This document is not intended for dissemination, publication or distribution to the public.  
It is provided for the exclusive use of licensed physicians.



B.1357.83.005.1

082400

*Creating Better Together™*